



DESCRIPTION AND INTENDED USE


MINI MITO BTE CLASSIC: Mini Mito, miniBTE digital hearing aid, battery 312. It is intended to compensate for mild to severe hearing losses. Easy to wear and comfortable.

Four (4) channels WDRC hearing aid, 12 adjusting bands (equalizers), 4 memories, feedback notch filter, low and high cut filters, output limiting. Moreover, Adaptive feedback cancelling and Noise reduction (10 bands and 3 levels of noise).

FEATURES

- 4-Channel Wide Dynamic Range Compression with Dynamic Contrast Detection™
- Patented Adaptive Feedback Cancelling
- Layered Noise Reduction™
- Expansion – Gain Reduction for Low-Level Sound
- Output Limiting – Wide Band Output Limiting (MPO) with Adjustable Thresholds
- Program Selector
- Battery Type: 312

FITTING RANGE

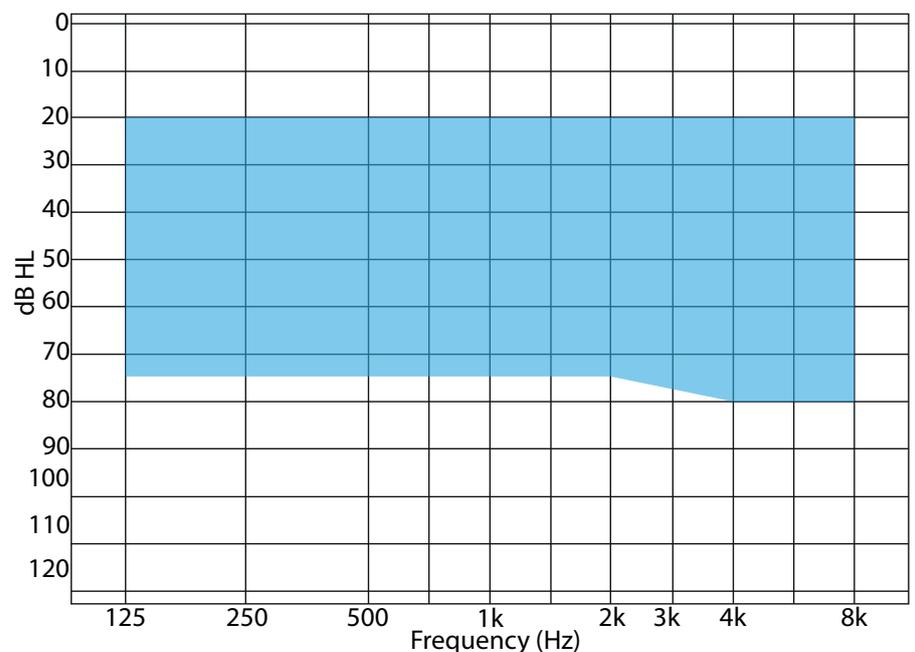

WARNING!
 This hearing aid can generate sound output levels in excess of 132 dB_{CP} (IEC 60318-4 Coupler).
 The hearing care specialist should be specially careful fitting the instrument as there may be risk of impairing the remaining hearing of the hearing instrument user.

CE FDA Approved
 1370

| Technical data <i>Measured according to</i> | | EAR SIMULATOR <i>IEC 60118-0:2015 and IEC 60318-5:2006</i> | 2CC Coupler <i>ANSI S3.22-2014, IEC60118-7:2005 and IEC 60318-5:2006</i> |
|---|---|--|--|
| Horentek Mini Mito | | CLASSIC | CLASSIC |
| Frequency range Hz | | 200/5500 Hz | 200/5500 Hz |
| OSPL90 | Peak Measurement | - | - |
| | Peak 1400 Hz | 126 dB | 121 dB |
| | HFA AVERAGE | 114 dB | 109 dB |
| Full-on gain 50* | Peak Measurement | - | - |
| | Peak 1400 Hz | 65 dB | 60 dB |
| | HFA AVERAGE | 43 dB | 38 dB |
| Reference test gain | | 45 dB | 35 dB |
| Telecoil output (1600 Hz) | 1 mA/m field 10 mA/m field SPLITS L/R | - | - |
| Total harmonic distortion (input 70 dB SPL) | 500 Hz | < 2 % | < 2 % |
| | 800 Hz | < 2 % | < 2 % |
| | 1600 Hz | < 2 % | < 2 % |
| Equivalent input noise level | | 31 dB SPL | 31 dB SPL |
| Battery consumption**(Battery 312) Typical | | 0.54 mA | 0.54 mA |
| Battery life, artificial measurements, hours *** | | 305 | 305 |

* Measured with the gain control of the hearing aid set to its full-on position minus 20dB and with an input SPL of 70 dB.
 ** Battery current is measured according to IEC 60118-0:2015 §7.7 after a settling time of a minimum of 3 minutes
 *** Based on the standardized battery consumption measurement. The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.

Data may vary by more or less 5%

